

Oncomine™ Reporter Digital Delivery and Installation

Pub. No. MAN0019656 Rev. E.0

This installation guide contains instructions for installing or updating the Oncomine™ Reporter Software using the web.

| | |
|--|---|
| ■ Before you begin | 1 |
| ■ Workflow for Digital Delivery | 1 |
| ■ Workflow for installation | 2 |
| ■ Download Oncomine™ Reporter Software | 2 |
| ■ Install Oncomine™ Reporter Software | 4 |
| ■ Activate license—Online | 4 |
| ■ Activate license—Offline | 5 |
| ■ Support | 7 |
| ■ Limited product warranty | 7 |

Before you begin

Ensure that you have command line access to the server where the Oncomine™ Reporter Software is being installed.

Workflow for Digital Delivery

Oncomine™ Reporter Software download workflow

You receive an email with link to sign in to the Publisher Licensing Center.

If you are signing in to the Publisher Licensing Center for the first time, update password.

Sign in to the Publisher Licensing Center.

Accept the terms and conditions.

Download Oncomine™ Reporter package.

Oncomine™ Reporter Software installation workflow


Copy and unpack Oncomine™ Reporter package.

Install Oncomine™ Reporter Software.

Activate software using the Activation ID.

Download Oncomine™ Reporter Software

1. Using the link from the welcome email, access the Publisher Licensing Center.

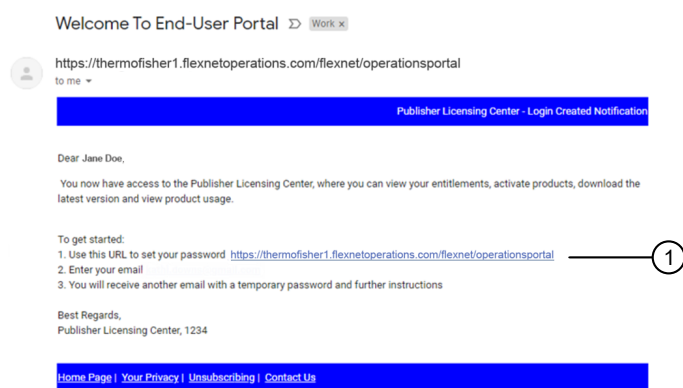


- ① Link to Publisher Licensing Center

- In the **Password Reset** screen, in the **User Id** box, enter your email, then click **Continue**.

The email should match the email in the welcome email.


A password reset email is sent to the email address that you entered.

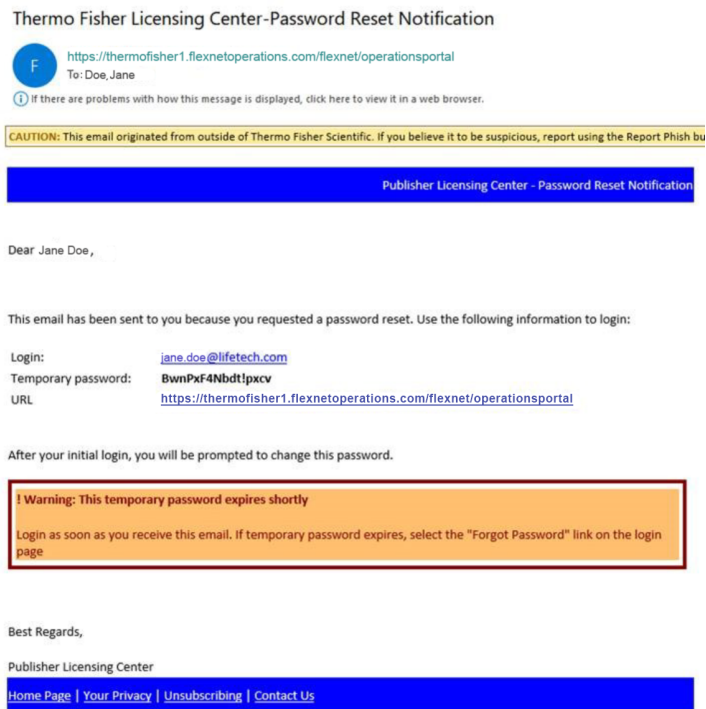

- Click the URL in the reset password email to open the **Change Password** screen.

- In the **Change Password** screen, in the **Password** box, enter the temporary password from the password reset email, enter a new password, confirm the password, then click **Continue**.

- Sign in to the Publisher Licensing Center using your credentials.

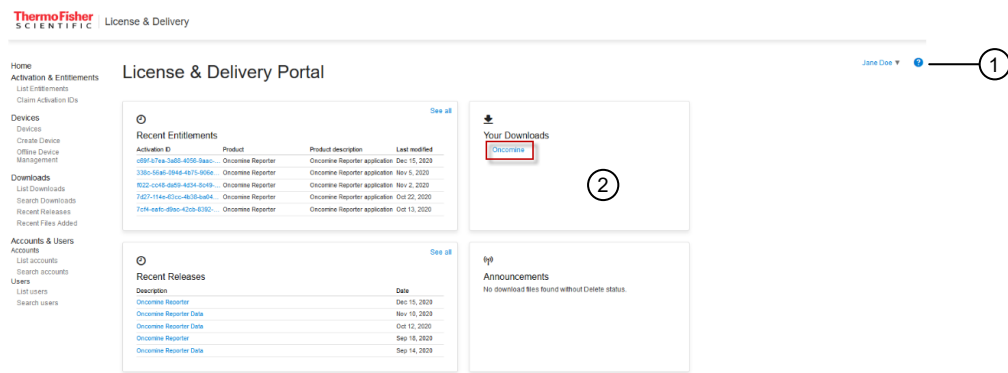

① Help button

② **Your Downloads** section

- (Optional) For help with the **License and Delivery Portal**, click ①.

- In the **License and Delivery Portal**, in the **Your Downloads** section, click **Oncomine**.

- In the **Download Packages** screen, in the **New Versions** tab, select **Oncomine Reporter**.

- In the **Software Terms and Conditions** screen, click **I agree**.

- In the **Downloads** screen, click on the file name to download the Oncomine™ Reporter ZIP file.

Install Oncomine™ Reporter Software

For additional details, see the *Oncomine™ Reporter User Guide* for your version of the Oncomine™ Reporter Software.

1. Copy the downloaded ZIP file to the server where you are installing the Oncomine™ Reporter Software.
2. Unpack the ZIP file into a folder. In the command screen, type the following: `unzip okr-5.3.zip`
3. Change to the OR5_3 directory that was created by unzipping the OKR 5.3.zip file: `cd ./OR5_3`
4. Run the installation script in the OR5_3 directory: `sudo bash ./install.sh`
5. After the new software is installed, select one of the following actions.
 - If the software is being updated, reload the page. The old URL becomes invalid and you must reload the page.
 - If you are installing the software for the first time, open a browser and navigate to the Oncomine™ Reporter.

Activate license—Online

1. Sign in to the Oncomine™ Reporter Software.
2. In the **Activate License** screen, enter the activation ID, then click **Activate**.
The activation ID is located in the confirmation email.

Entitlement details

Dear Jane Doe,

Thank you for your recent purchase. The details listed below contain information to help you get started using your products and services which are immediately available. For software purchases, please follow the instructions at the end of this email to begin using your product(s).

| Order Details | | | |
|----------------|---|-------------|--------------|
| Order Number: | 5f14-c833-4bfa-4850-9457-dd89-57ea-2c6d | Start Date: | Dec 15, 2020 |
| Product Name: | Oncomine Reporter | Exp. Date: | Dec 15, 2021 |
| Activation ID: | c69f-b7ea-3a88-4056-9aac-54b4-39c7-cf9c | Quantity: | 1 |

Product Download and Licensing Information

All related software is available for immediate download via the [Download and License Center](#) available 24/7. The email address associated with the order will be required to access the site.


This site also contains all relevant documentation and license activation codes (if applicable) for your recent purchase.

Customer Service Information

Along with the purchase of your software you also have full access to our customer support portal. Here you will be able to find FAQ, Knowledge Base as well as gain access to the [Customer Community](#)

If there are any questions in regards to your order please [contact us](#) for additional assistance

[Home Page](#) | [Your Privacy](#) | [Unsubscribe](#) | [Contact Us](#)

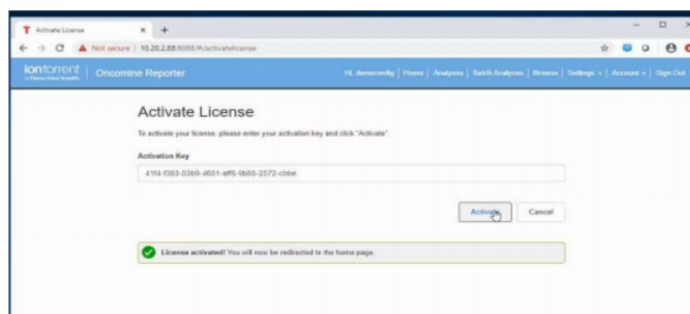

The software is ready to use.

Activate license—Offline

1. Sign in to the Oncomine™ Reporter Software.
2. In the **Activate License** screen, enter the activation ID, then click **Download**.

The activation ID is located in the confirmation email.

Entitlement details

Dear Jane Doe,

Thank you for your recent purchase. The details listed below contain information to help you get started using your products and services which are immediately available. For software purchases, please follow the instructions at the end of this email to begin using your product(s).

| Order Details | | | |
|----------------|---|-------------|--------------|
| Order Number: | 5f14-c833-4bfa-4850-9457-dd89-57ea-2c6d | Start Date: | Dec 15, 2020 |
| Product Name: | OncoPrint Reporter | Exp. Date: | Dec 15, 2021 |
| Activation ID: | c69f-b7ea-3a88-4056-9aac-54b4-39c7-cf9c | Quantity: | 1 |

Product Download and Licensing Information

All related software is available for immediate download via the [Download and License Center](#) available 24/7. The email address associated with the order will be required to access the site.

This site also contains all relevant documentation and license activation codes (if applicable) for your recent purchase.

Customer Service Information

Along with the purchase of your software you also have full access to our customer support portal. Here you will be able to find FAQ, Knowledge Base as well as gain access to the [Customer Community](#)

If there are any questions in regards to your order please [contact us](#) for additional assistance

[Home Page](#) | [Your Privacy](#) | [Unsubscribe](#) | [Contact Us](#)

Activate License

OncoPrint Reporter wasn't able to contact the license server. To activate your license, please follow these steps:

Step 1: Enter your activation ID and click "Download" to download an activation file (okr-activation.bin).

Activation ID

XXXX-XXXX-XXXX-XXXX-XXXX-XXXX-XXXX

Download

Step 2: Go to <https://thermofisher1.flexnetoperations.com/flexnet/operationsportal>, click on "Offline Device Management", upload the activation file you downloaded in step 1 (okr-activation.bin), and download the license file (capabilityResponse.bin).

Step 3: Return to this page. Select the license file downloaded in step 2 (capabilityResponse.bin) and click "Activate".

Upload license file




Activate

Cancel

The offline activation file is generated and saved to your local storage (okr-activation.bin).

- Open a new browser window and access the **License & Delivery Portal** (<https://thermofisher1.flexnetoperations.com/flexnet/operationsportal>).

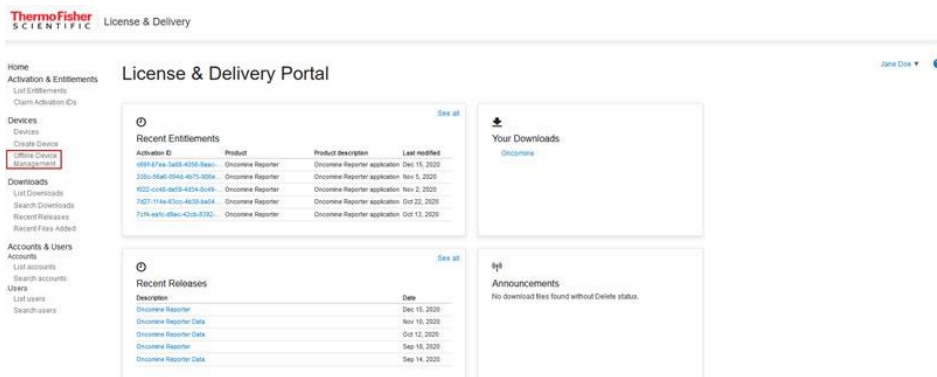

- In the left navigation pane, under the **Devices** section, click **Offline Device Management**.
- In the **Offline Device Management** screen, ensure that **Generate license** is selected, click **Choose File**, then select the activation file downloaded in step 2 (okr-activation.bin).
- Click **Upload**.
The license file is generated and saved to your local storage (capabilityResponse.bin).
- In the browser window with the **Activate License** screen, in step 3 in the **Activate License** screen, upload the license file from step step 6, then click **Activate**.

Support

For customer and technical support, contact your local support representative.

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale at www.thermofisher.com/us/en/home/global/terms-and-conditions.html. If you have any questions, please contact Life Technologies at www.thermofisher.com/support.



Life Technologies Corporation | 5781 Van Allen Way | Carlsbad, California 92008 USA
For descriptions of symbols on product labels or product documents, go to thermofisher.com/symbols-definition.

The information in this guide is subject to change without notice.

DISCLAIMER: TO THE EXTENT ALLOWED BY LAW, THERMO FISHER SCIENTIFIC INC. AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

Revision history: Pub. No. MAN0019656

| Revision | Date | Description |
|----------|-------------------|---|
| E.0 | 13 September 2021 | Updated to correspond to Oncomine™ Reporter 5.3. |
| D.0 | 11 June 2021 | Updated to correspond to Oncomine™ Reporter 5.2. |
| C.0 | 15 April 2021 | Updated to correspond to current process for the installation and activation of Oncomine™ Reporter Software using digital delivery. |
| B.0 | 4 March 2021 | Updated to correspond to Oncomine™ Reporter 5.1. |
| A.0 | 17 December 2020 | New document—Installation and activation of Oncomine™ Reporter Software using digital delivery. |

NOTICE TO PURCHASER: DISCLAIMER OF LICENSE: Purchase of this software product alone does not imply any license under any process, instrument or other apparatus, system, composition, reagent or kit rights under patent claims owned or otherwise controlled by Thermo Fisher Scientific, either expressly, or by estoppel.

©2021 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.